



# NEW Model NU-99100J



Ultralow Under-Counter Freezer

-86°C ULTRALOW TEMPERATURE FREEZER

The NuAire Blizzard -86°C Ultralow Mini-Freezer features reliable performance for the storage of stem cells, blood & tissue samples, and vaccine vials while occupying the smallest space among our ultralow freezers. The NU-99100J model can maintain temperature uniformity of  $\pm 5^{\circ}\text{C}$  at  $-80^{\circ}\text{C}$  and features a simplified LED display screen, a quadruple gasketed outer door, 2 interior steel mini-doors, and produces visual/audible alarms. The quiet compressor motor utilizes non-toxic Hydrocarbon refrigerants and operates at a quiet 47 dB(A) while maximizing energy efficiency.



## Storage Capacity: 100 L (3.5 ft<sup>3</sup>)

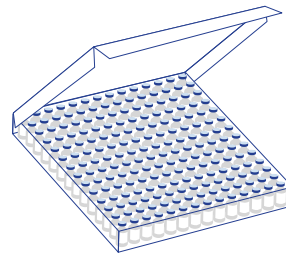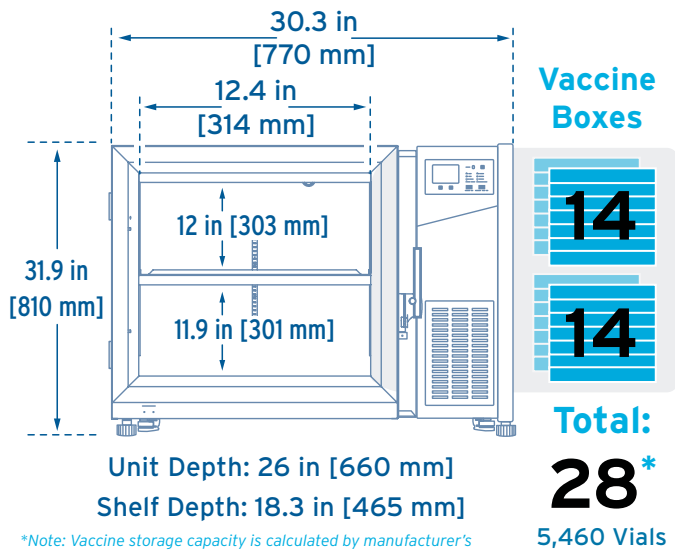

1 Box = 195 Vials

Box Width: 9" [229 mm]  
Box Depth: 9" [229 mm]  
Box Height: 1.6" [41 mm]



1 Vial = 5 Doses

\*Note: 2 Doses are needed to inoculate 1 person

Model	Configuration	Storage Capacity	Interior Dimensions (W x D x H)	Exterior Dimensions (W x D x H)	Recommended Temperature Range	Net Weight	Electrical Configurations
NU-99100J	Under-Counter	3.5 cu ft [100 L]	13 x 18.9 x 24.8 in [330 x 481 x 630 mm]	30.2 x 26 x 31.9 in [770 x 660 x 810 mm]	-50°C to -86°C	238 lbs [108 kg]	115 VAC/60 Hz 230 VAC/50 Hz



Click or Scan QR Code to Learn More



NuAire, Inc. | 2100 Fernbrook Lane | Plymouth, MN 55447 | U.S.A. | [www.nuaire.com](http://www.nuaire.com)

20-1583-F-G-EN-2-0121 © Copyright 2021. NuAire, Inc. All Rights Reserved.